waterpik™

Waterpik™ Sensonic™ Toothbrush STW Series

Waterpik™ Sensonic™
Zahnbürste der STW-Serie

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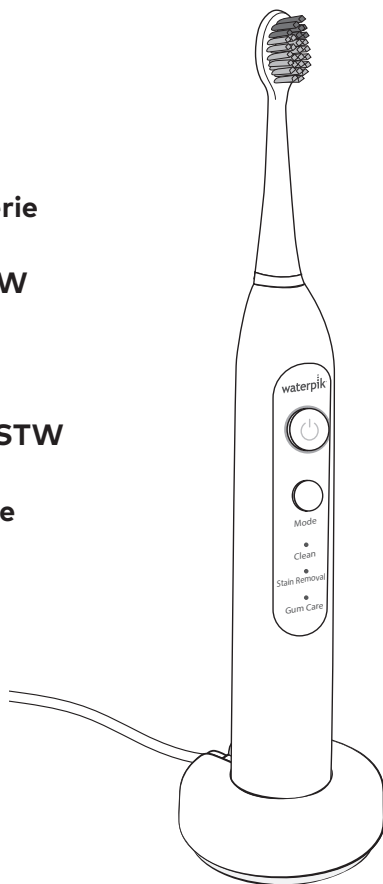
SCALE: 100%
UPC: N/A
SIZE: 10" x 7" flat; 5" x 7" folded
FORM NUMBER: 20033038-F AA
INK: Process Black
NOTES:

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not handle the user supplied USB power supply (power adapter) with wet hands.
- Unplug the user supplied USB power supply (power adapter) when not in use.
- Check the charging cord and user supplied USB power supply (power adapter) for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only for its intended use as described in these instructions or as recommended by your dental professional.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact your local distributor which can be found at www.waterpik.com/intl
- A user supplied USB-A wall charger (power adapter), safety certified (IPX4, Class II) with a rated output of 5.0 VDC, 1.0 A Max is required to charge this product.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Keep the charging cord away from heated surfaces.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

- Only use replacement brush heads and accessories that are recommended by Water Pik, Inc. And as identified in this instruction manual.
- If the charging cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- There are no consumer serviceable electrical parts in this product and it does not require maintenance.
- Ensure contacts on the cord, the unit and the power adapter are clean and dry while the unit is being charged.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet. Contact your doctor with any questions.
- Do not expose the battery or appliance to excessive temperatures.
- Using toothpastes containing peroxide may damage the plastic housing of this product.

SAVE THESE INSTRUCTIONS

TABLE OF CONTENTS



Product Description	5
Getting Started with Your Sensonic™ Toothbrush	6
Care of Your Sensonic™ Toothbrush	7
Battery Care, Battery Disposal/Recycling and Battery Removal	8
Warranty Information	10

SYMBOLS			
The following symbols may appear on your product and/or packaging. Some of these symbols may not be relevant in your region and are listed for information purposes only.			
Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Do not immerse this device in any liquid.
	Manufacturer		Do not use this device in a shower, bathtub and do not immerse in any liquid.
IPX7 or IPX4	Ingress Protection Code		Authorized Representative in the European Community
	IEC Class II Equipment		CE marking of conformity
	Detachable Power Supply		UK Conformity Assessed
	Direct Current (DC) Voltage		Regulatory Compliance Mark
	Alternating Current (AC) Voltage		Separate collection for waste electrical and electronic equipment.


PRODUCT DESCRIPTION

Sensonic™ Toothbrush

1. Sensonic™ Toothbrush Brush Head
2. Toothbrush power button
3. LED charge indicator
4. Toothbrush mode button
5. LED mode indicator
6. Charging puck
7. USB charging cable
8. Toothbrush travel case

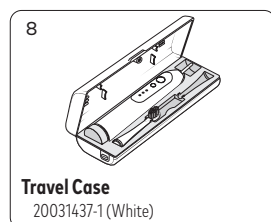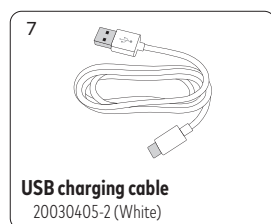
**Brush Heads****Sensonic™ Toothbrush Brush Heads***

20030861-3

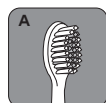
*Replace Every 3 months

Replacement brush heads/accessories can be purchased from your local distributor, which can be found at www.waterpik.com/intl.

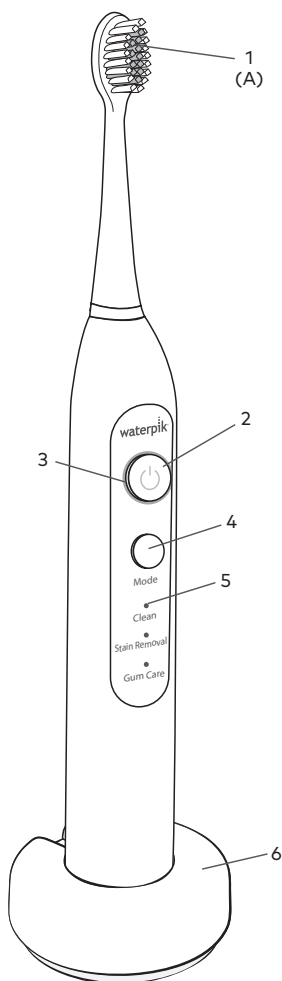**Sensonic™ Toothbrush**

GETTING STARTED WITH YOUR SENSONIC™ TOOTHBRUSH

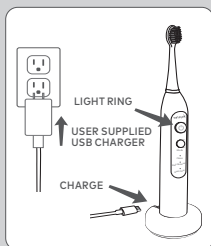
Charging the Sensonic™ Toothbrush

Charge for 4 hours before first use. Light ring will flash on handle while charging.

When fully charged all lights will shut off. Light ring will glow red when battery is running low. A full charge will last for approximately 2 weeks of typical use (2 minutes of brushing twice a day).

Note: Refer to *IMPORTANT SAFEGUARDS* section for required USB wall charger specifications.

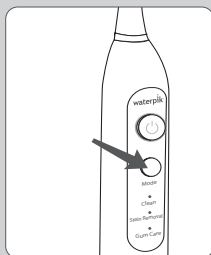Modes



Clean: Outstanding everyday cleaning performance.

Stain Removal: Gently removes surface stains.

Gum Care: Gently enhances gum stimulation.

Sensonic™ brushing

- Slide the brush head onto handle and push it down firmly.
- Place the brush head into your mouth before turning the unit on.
- Select mode by pressing mode button.
- For best results, position brush along gumline at a 45-degree angle.
- Use light pressure with a small, slow back and forth motion covering one or two teeth at a time.
- **Quad Timer:** The Toothbrush will create a short pause every 30 seconds, to remind you to move to the next quadrant of your mouth.
- **2-Minute Timer:** The brush will automatically shut off after 2 minutes of run time.

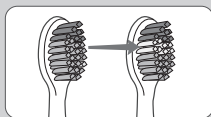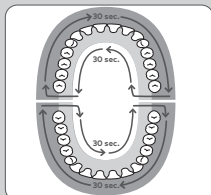Bristles change from green to yellow when time to replace. (For optimal performance replace brush heads every 3 months.)

CARE OF YOUR SENSONIC™ TOOTHBRUSH

**Cleaning**

1. **After each use**, turn off unit, remove the brush head from the handle and rinse brush head under running water.
2. Wipe the handle dry after each use.
3. Periodically, clean the charger puck, first unplugging the unit and then wiping with a damp cloth.

Battery Care

If the toothbrush is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage.

Do not expose the battery or appliance to excessive temperatures.

Battery Disposal/Recycling





This product contains a **rechargeable lithium-ion battery**. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations.

Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

Please contact your local waste administrative department for any questions.

NOTE: The rechargeable battery is not replaceable, remove only at end of product life. This process is not reversible and voids the product warranty.

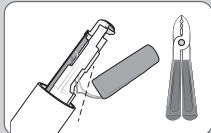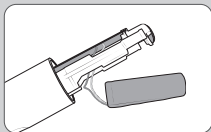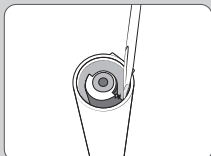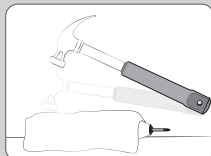Battery Removal

To remove the rechargeable battery for recycling, please discharge the battery fully before removing.

Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work.

1. Run down the battery by operating the BRUSH button until the motor does not run.
2. Obtain a towel or cloth, hammer, flat-head screwdriver, and wire cutters.
3. Remove and discard the brush head. Cover the handle entirely with the towel or cloth.
4. Clamp the top of the handle against your work surface with one hand and strike the handle housing around the bottom perimeter on all sides to release the end cap latch. Continue striking and rotating the handle until the end cap releases.
5. Remove the internal components by turning the handle upside down, then un-latch the plastic tabs and slide the chassis from the outer housing.
6. The internal components should slide out when the handle is positioned upside-down. If this does not occur, repeat steps 5 and 6 until the components release.
7. Remove the battery from the chassis. Cut the battery wires one at a time with the wire cutters.
8. Recycle and dispose of the battery properly. The remaining product pieces should be discarded appropriately.

WARNING:

DO NOT cut both battery wires at the same time.

DO NOT disassemble, puncture, hit, squeeze or throw the battery into fire.

WARRANTY INFORMATION

Service Maintenance

Waterpik™ Sonic™ Toothbrushes have no consumer - serviceable electrical items and do not require routine service maintenance. For all your service needs go to www.waterpik.com/intl to find your local distributor. Some replacement parts and accessories for the toothbrush are also available from your local distributor which can be found at www.waterpik.com/intl.

Refer to the serial and model numbers in all correspondence. These numbers are located on the base or handle of your toothbrush.

Still have questions? Go to www.waterpik.com/intl to find your local distributor.

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase (Spain Only: 3 years). Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips, etc. **Use with improper voltage will damage product and VOID warranty.**

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

To locate your local distributor, go to www.waterpik.com/intl.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

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©2024 Water Pik, Inc.
**Water Pik, Inc., a subsidiary
of Church & Dwight Co., Inc.**
SOFIBEL SAS