Waterpik® Water Flosser Models WF-11 | WF-12

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
Idropulsore Waterpik® Modelli WF-11 | WF-12

Irrigador bucal Waterpik®
Modelos WF-11 | WF-12

Waterpik® Water Flosser
Modellen WF-11 | WF-12

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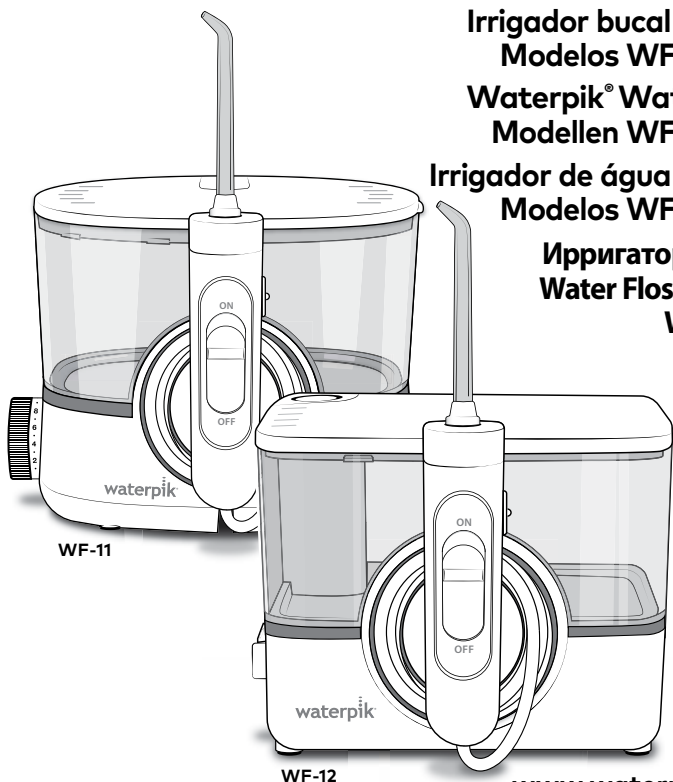DO NOT PRINT

DATE: 11.21.2022
CLIENT: WATERPIK

FILENAME: 20029228-FAB_WF-11-12_EU_IM_F.indd
DESCRIPTION: WF-11-12 Instruction Manual
LANGUAGES: 8 (Eng, Ger, Fr, It, Span, Dut, Port, Russ)
FILE TYPE: Adobe InDesign CC
SCD JOB# / PREPARED BY: 22WPOH7491 / BT / BB
CLIENT ROUND: F

SCALE: 100%
UPC: N/A
SIZE: 10" x 7" flat; 5" x 7" folded
FORM NUMBER: 20029228-F AB
INK: Process Black
NOTES:

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not handle the user supplied USB charger with wet hands.
- Do not place or drop the user supplied USB charger into water or other liquid.
- Do not use this product while bathing.
- Do not place or store this product where it can fall or be pulled into a tub or sink.
- Do not reach for this product if it has fallen into water. Unplug from electrical outlet immediately.
- Check the charging cord and user supplied USB charger for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Use this product only as indicated in these instructions or as recommended by your dental professional.
- A user supplied USB charger with the following specifications is required to charge this product: USB Type-A, Rated Output: 5.0 VDC, 1.0 A.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Only use tips and accessories recommended by Water Pik, Inc.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet.
- Do not drop or insert any foreign object into any opening or hose.
- Keep power supply/charger cord away from heated surfaces.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

- Do not use iodine or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product. Use of these products may void warranty.
- Fill reservoir with water or other dental professional recommended solution only.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip (PP-100E).
- Do not use this product for more than five minutes in each two-hour period.
- Unit will shut off automatically after three minutes of continuous usage.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the supply cord is damaged, it must be replaced by a special cord set available from the manufacturer or its service agent.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. at 1-800-525-2774 or visit www.waterpik.com/intl.

SAVE THESE INSTRUCTIONS


TABLE OF CONTENTS

Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik® Water Flosser	8
Battery Removal	9
Warranty Information	10

SYMBOLS			
The following symbols may appear on your product and/or packaging. Some of these symbols may not be relevant in your region and are listed for information purposes only.			
Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Do not immerse this device in any liquid.
	Manufacturer		Do not use this device in a shower, bathtub and do not immerse in any liquid.
IPX7 or IPX4	Ingress Protection Code		Authorized Representative in the European Community
	IEC Class II Equipment		CE marking of conformity
	Detachable Power Supply		UK Conformity Assessed
	Direct Current (DC) Voltage		Regulatory Compliance Mark
	Alternating Current (AC) Voltage		Separate collection for waste electrical and electronic equipment.

PRODUCT DESCRIPTION

1. Power ON/OFF
2. Handle
3. Water ON/OFF switch
4. Pressure control
5. Magnetic handle cradle
6. Handle Hose Swivel
7. Tip
8. Tip rotator
9. Reservoir lid
10. Reservoir
11. Tip eject button
12. Charging Contacts
13. Base

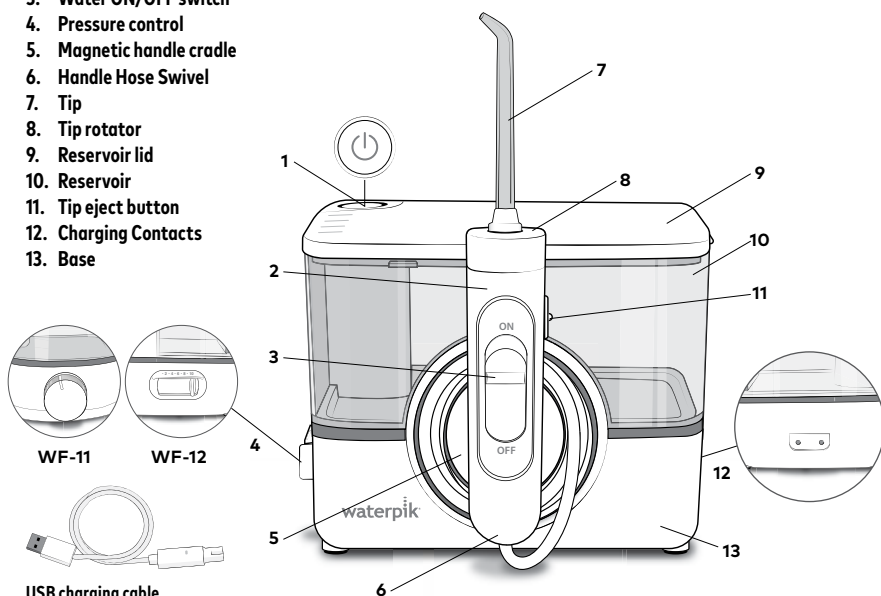
USB charging cable
20030698-1 White
20030698-2 Black

TIPS






- A. Orthodontic Tip***
- Braces/General Use
- B. Pik Pocket™ Tip***
- Periodontal Pockets/Furcations
- C. Plaque Seeker™ Tip***
- Implants/Crowns/Bridges/
Retainers/General Use
- D. Classic Jet Tip****
- General Use
- E. Implant Denture Tip****
- Fixed Implant Bridges and Dentures
- F. Tongue Cleaner****
- Fresher Breath

NOT ALL TIPS INCLUDED IN ALL MODELS

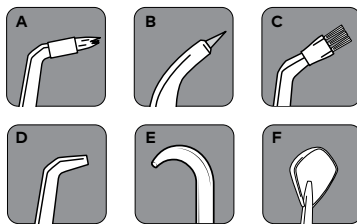*Replace Every 3 months **Replace Every 6 months

Replacement tips can be purchased from your local distributor, which can be found at www.waterpik.com/intl.

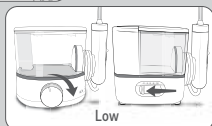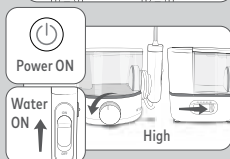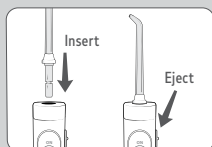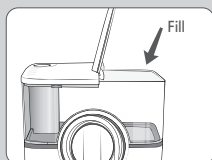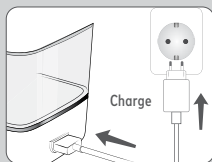Charging Unit Prior to First Use

A user supplied USB charger with the following specifications is required to charge this product: USB Type-A, Rated Output: 5.0 VDC 1.0 A. **Ensure contacts on BOTH the cable and unit are clean and dry before connecting the charger.** Plug the USB-charger into a wall outlet and connect the magnetic charger cable to the unit. The charge indicator on the unit will glow while charging and then begin to flash momentarily every 3 seconds when the unit is fully charged. A full charge will last for approximately 4 weeks when unit is used once a day for one minute a day. When the batteries are in need of charging (approximately 30% charge left) the charge indicator light will blink 3 times quickly to indicate that the user should consider recharging the unit. A complete recharge will take 4-5.5 hours. If you use your Water Flosser once a day or less, it is recommended that you charge it every 3-4 weeks. If you use it twice a day or more, charge the unit more frequently.

Preparing the Reservoir

Remove the reservoir, lift the reservoir lid, and fill with warm water. Fill with warm water. Replace the reservoir over the base and press down firmly and close the lid.

Inserting and Removing Tips

Insert tip into the center of the Water Flosser handle, pressing down firmly until the tip clicks into place. The colored ring will be flush with the top of the handle if the tip is correctly locked into place. To remove tip from the handle slide the tip eject button on the handle.

Prime Unit (First Time Use Only)

Set pressure control on base to high (pressure setting No.10), point the handle and tip into sink. Turn power and water ON until water flows from the tip.



Adjusting the Pressure Setting

Set the pressure control on the base of the unit to the lowest setting (pressure setting No.1) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

Using the Water Flosser

Always place the tip in your mouth before turning the unit power ON. With water switch on handle in OFF position, turn power on. **Lean over sink and slightly close lips enough to prevent splashing while still allowing water to flow from mouth into the sink. Turn water switch ON. When finished, turn water switch and power OFF.** During use you can stop the flow of water by sliding the water switch on handle to the OFF position.

NOTE: Motor will continue to run.

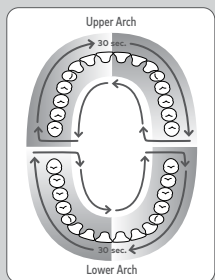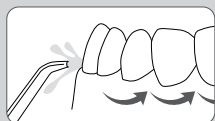For best results, start with back teeth and work toward the front teeth. Direct the stream at a 90-degree angle to your gumline with tip in mouth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

One-Minute Timer/30-second Pacer

The Water Flosser will create a short pause at 30-seconds, to remind you to move to the other arch of your mouth and another short pause at 1-minute to signal you have reached the recommended cleaning time. The unit will continue to run without pauses until the unit is turned off. Depending on the pressure setting, some water will remain in the reservoir after 1-minute.

When Finished

Turn the power OFF. Empty any liquid left in the reservoir.

NOTE: If you used mouthwash or other solution besides water, see "Using Mouthwash and Other Solutions" section.

NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.


To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting.** Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.




Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

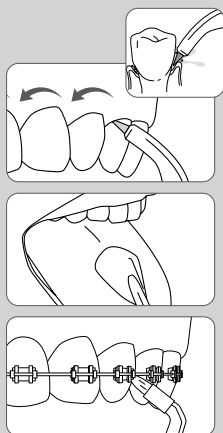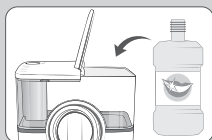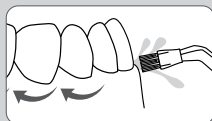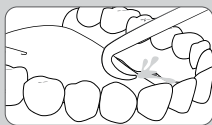Implant Denture Tip

The Implant Denture tip is specifically designed to access areas that are unreachable with a Classic Jet tip, Plaque Seeker™ tip, Orthodontic tip, or Pik Pocket™ tip. To use the Implant Denture tip, place the tip on the lingual side at one end of the implant retained denture. The opening of the tip is directed between the natural gum tissue and the denture directing the water under the denture. Slowly glide the tip along the space between the natural gum tissue and the denture until you have cleaned underneath the entire length of the denture.

Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Before cleaning, disconnect power supply/charger cable. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

The water reservoir is removable for easy cleaning. Place the reservoir in the top rack of your dishwasher for cleaning. We recommend air drying instead of a heated dry cycle.


Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 2 tablespoons of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.

Leakage

If leaking due to a damaged hose, replacement hoses and handles can be purchased from your local distributor, which can be found at www.waterpik.com/intl.

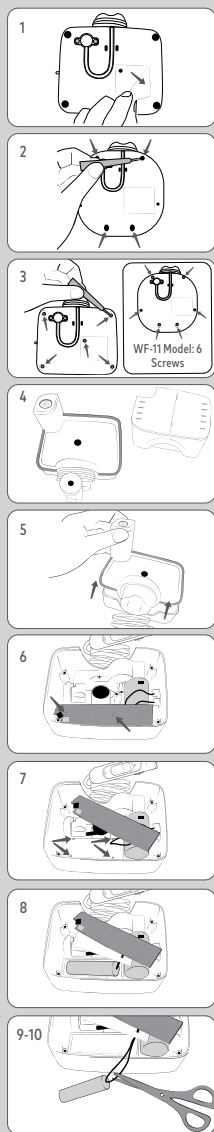Battery Removal

The rechargeable lithium-ion battery inside this product cannot be replaced and shall be removed at the end of the product's life prior to recycling.

This process is not reversible and voids the product warranty.

Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and the surface on which you work. Deplete the battery charge until the unit no longer operates prior to battery removal.

1. Turn unit over. Remove product label (Model WF-12 Only)
2. Remove the 4 rubber feet by carefully prying them upward with a small flat head screwdriver.
3. Remove 5 screws from the indicated areas with a small Phillips head screwdriver.
4. Turn unit upright and remove reservoir and reservoir lid.
5. Remove upper housing by pulling it upward as shown.
6. Remove screws from the indicated areas with a small Phillips head screwdriver.
7. Pull circuit board assembly aside. Remove 4 screws from the indicated areas with a small Phillips head screwdriver.
8. Lift up battery cover as shown.
9. Remove battery from housing.
10. Cut the black battery wire with a pair of scissors or side-cutters then cut the red battery wire.

NOTE: DO NOT CUT BOTH BATTERY WIRES AT THE SAME TIME.

11. Please recycle the unit and safely dispose of the battery per local or state requirements.

Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your warranty and accessory needs please contact your local distributor. Find your local distributor at www.waterpik.com/intl. Please do not ship product to Water Pik, Inc. offices, this will delay service. Refer to the serial and model numbers in all correspondence. These numbers are located on the back of the handle and bottom of the reservoir respectively.

Still have questions? Go to www.waterpik.com/intl to find your local distributor.

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years (Spain Only: 3 years), from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Repair or replacement services provided for by the guarantee will be carried out by the local distributor, whose details may be found on www.waterpik.com/intl. Upon request of the consumer, a copy of the guarantee complete with the local distributor's information will be provided by email or mail. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary depending on the law in your place of residence. To locate your local distributor, go to www.waterpik.com/intl.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

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©2022 Water Pik, Inc.
**Water Pik, Inc., a subsidiary
of Church & Dwight Co., Inc.**
Form No. 20029228-F AB