waterpik™

Waterpik™ Water Flosser Model WF-600 Series

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DO NOT PRINT

DATE: 10.28.24
CLIENT: WATERPIK

FILENAME: 20033681-FAB_WF-600_Series_EU_IM_F.indd
DESCRIPTION: WF-600 Series Instruction Manual
LANGUAGES: 7 (Eng, Ger, Fr, It, Span, Dut, Russ)
FILE TYPE: Adobe InDesign CC
SCD JOB# / PREPARED BY: 24WPOH8389 / BT
CLIENT ROUND: F

SCALE: 100%
UPC: N/A
SIZE: 10" x 7" flat; 5" x 7" folded
FORM NUMBER: 20033681-F AB
INK: Process Black
NOTES:

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop the product or external power supply into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not handle the external power supply with wet hands.
- Check the supply cord and external power supply for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only for its intended use as described in this Instruction Manual or as recommended by your dental professional.
- Never operate this product if it has a damaged cord or external power supply, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact your local distributor which can be found at www.waterpik.com/intl.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or external power supply.
- Keep electrical cord away from heated surfaces.
- When unplugging the external power supply, be sure to pull the power supply itself. **Do not** pull the cord to unplug.
- Never drop or insert any object into any hose or opening.

- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Remove water flosser reservoir to fill. Fill reservoir with warm water or other dental professional recommended solutions only.
- Only use tips and accessories recommended by Water Pik, Inc. and as identified in this instruction manual.
- The supply cord and external power supply cannot be replaced. If the cord or external power supply are damaged, the appliance must be scrapped.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Do not direct water from the water flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist and/or medical professional before using the water flosser or any other oral hygiene aid.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- Do not use this product for more than five minutes in each two-hour period.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
- This product does not contain any user serviceable electrical parts and it does not require maintenance. Do not attempt to open the product.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device containing magnets or electronics, contact your physician or the device manufacturer prior to using the appliance.

EU EcoDesign power per (EU) 2023/826

- Operating Modes = Flossing and Standby
- Standby Mode power consumption = 0.4 W
- Time to Standby mode; User shutoff = 0 minutes
- Time to Standby mode; Auto-shutoff = 5 minutes

SAVE THESE INSTRUCTIONS


TABLE OF CONTENTS

Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik™ Water Flosser	9
Warranty Information	10

SYMBOLS GLOSSARY (Not all symbols apply to all products)			
Explanation of symbols used on Water Pik, Inc. product labels and packaging. Applicable symbols for a specific product are found on that product's labels and packaging.			
Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Do not immerse this device in any liquid.
	Manufacturer		Do not use this device in a shower, bathtub and do not immerse in any liquid.
IPX7 or IPX4	Ingress Protection Code		Authorized Representative in the European Community
	IEC Class II Equipment		CE marking of conformity
	Detachable Power Supply		UK Conformity Assessed
	Direct Current (DC) Voltage		Regulatory Compliance Mark
	Alternating Current (AC) Voltage		Separate collection for waste electrical and electronic equipment.
	Contains a recyclable Li-ion battery		Contains a recyclable Ni-MH battery

PRODUCT DESCRIPTION


1. Dual purpose covered lid
2. Reservoir
3. Tip storage under lid
4. Power button
5. Mode selector button
6. LED mode indicators
7. Pressure control dial
8. Exchangeable tips
9. Rotating control knob
10. Tip eject button
11. Water on/off switch
12. Handle
13. Water hose
14. Base
15. External power supply

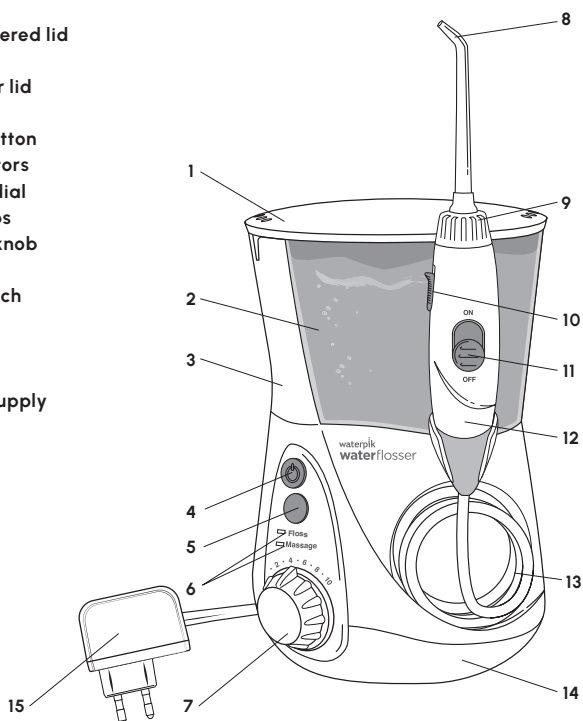TIPS

- A. **Orthodontic Tip***
- Braces/General use
- B. **Pik Pocket™ Tip***
- Periodontal pockets/
Furcations
- C. **Plaque Seeker™ Tip***
- Implants/Crowns/Bridges/
Retainers/General use
- D. **Classic Jet Tip****
- General use
- E. **Tongue Cleaner****
- Fresher breath





NOT ALL TIPS INCLUDED IN ALL MODELS

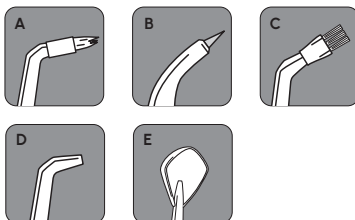*Replace every 3 months **Replace every 6 months

Replacement tips can be purchased from your local distributor, which can be found at www.waterpik.com/intl.

GETTING STARTED

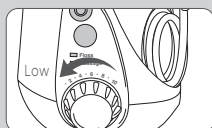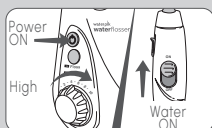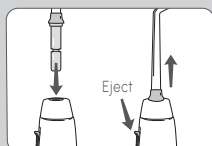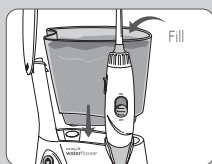**PREPARING THE UNIT****Powering the Unit**

Plug the external power supply into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

Preparing the Reservoir

Lift the reservoir lid, remove the reservoir, and fill with lukewarm water. Replace the reservoir over the base and press down firmly and close the lid.

Inserting and Removing Tips

Insert tip into the center of the water flosser handle, pressing down firmly until the tip clicks into place. The colored ring will be flush with the top of the handle if the tip is correctly locked into place.

To remove tip from the handle push the tip eject button on the handle.

Caution: Do not eject tip while unit is running.

Prime Unit (First Time Use Only)

Set pressure control dial on base to high (pressure setting No.10), point the handle and tip into sink. Turn power ON and water ON until water flows from the tip. Turn water OFF and power OFF.

Adjusting the Pressure Setting

Turn the pressure control dial on the base of the unit to the lowest setting for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

Using the Water Flosser

Always place the tip in your mouth before turning the water ON. With water switch on handle in the OFF position, turn power ON.



Lean over sink and slightly close lips enough to prevent splashing while still allowing water to flow from mouth into the sink. Turn water switch ON. When finished, turn water switch and power OFF. During use you can stop the flow of water by sliding the water switch on handle to the OFF position.

NOTE: Motor will continue to run.

Recommended Technique

Lean low over sink and place tip in mouth. Direct the jet stream at a 90-degree angle to your gumline. For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth. Slightly close lips to avoid splashing but allow water to flow freely from mouth into the sink.

Tip: LOOK INTO SINK, not in mirror, to reduce splash.

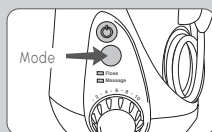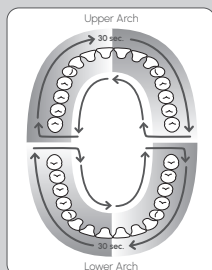One-Minute Timer/30-second Pacer (Floss Mode Only)

The water flosser will create a short pause at 30-seconds, to remind you to move to the other arch of your mouth and another short pause at 1-minute to signal you have reached the recommended cleaning time. The unit will continue to run without pauses until the unit is turned off. Depending on the pressure setting, some water will remain in the reservoir after 1-minute.

Timer may be deactivated if desired. Steps to deactivate:

1. While unit is off, press the mode button once to wake the unit up. One of the LEDs should light up.
2. Press and hold the mode button for approximately 3-seconds until the **Massage Mode** light blinks. The timer is now deactivated.

If you wish to reactivate the timer:

1. While unit is off, press the mode button once to wake the unit up. One of the LEDs should light up.
2. Press and hold the mode button for approximately 3-seconds until the **Floss Mode** light blinks. The timer is now activated.

Or unplug unit to reset to factory settings.

NOTE: If one of the LED's is on and you hold the mode button down for 3 seconds, one of the LED's will blink and the stutter feature will toggle.

Hydro-Pulse Massage Mode




Hydro-Pulse Massage offers a unique water pulsation pattern designed to enhance gum stimulation. For best results, use the Floss Mode as instructed. After using Floss Mode for 1-minute, switch the unit to Massage Mode, adjust pressure as desired, and aim the water stream at your gums and massage throughout your mouth.

To determine in which mode the unit is set while it is off simply press mode button and the LED light will indicate current mode. To change mode, press mode button again. You may also press the mode button while unit is on to change mode during use.

When Finished

Turn the water OFF and unit OFF. Empty any liquid left in the reservoir. Unplug the unit.

NOTE: If you used a mouthwash or other solution besides water, see *Using Mouthwash or Other Solutions* section.

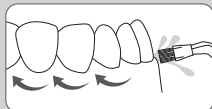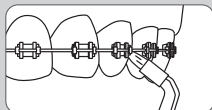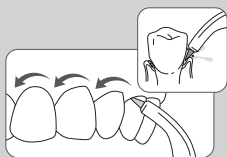NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

SPECIALTY TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets. To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and continue tracing along the gumline.

Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

Tongue Cleaner Tip

To use the tongue cleaner tip, set to lowest pressure setting, place tip in the center/middle of your tongue about half way back and start the water flosser. Pull forward with light pressure. Increase pressure as you prefer.


Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

Accessory Storage

Accessory tip storage for your two most frequently used tips is under the reservoir lid and can be easily accessed by removing the reservoir.

CLEANING AND TROUBLESHOOTING YOUR WATERPIK™ WATER FLOSSER

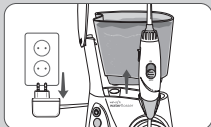Using Mouthwash and Other Solutions

Your Waterpik™ water flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior. The water reservoir is removable for easy cleaning. It is **top rack dishwasher safe with heat cycle OFF**. It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 part white vinegar and 2 parts warm water inside reservoir. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water. This process should be done every **1 to 3 months** to ensure optimal performance. If using mouthwash or other solutions, more frequent cleaning is recommended.

Leakage

If leaking due to a damaged hose, replacement hoses and handles are available from your local distributor which can be found at www.waterpik.com/intl.

NOTE: Ensure you are not storing any leftover water in your reservoir - this may cause leaking.

WARRANTY INFORMATION

Service Maintenance

Waterpik™ water flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all parts and accessories go to www.waterpik.com/intl. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit. Still have questions? Go to www.waterpik.com/intl.

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase (Spain Only: 3 years). **Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim.** Water Pik, Inc. will replace the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips, etc. **Use with improper voltage will damage product and VOID warranty.** This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

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©2025 Water Pik, Inc.
**Water Pik, Inc., a subsidiary of
Church & Dwight Co., Inc.**
SOFIBEL SAS
Form No. 20033681-F AB